

## MODEL

Jazzy Air 2	Jazzy Air 21001
Weight Capacity	300 lbs.
Turning Radius	24.25"
Overall Width <sup>2</sup>	25.25"
Overall Length <sup>2</sup>	42"
Maximum Speed <sup>1</sup>	Up to 4 mph (elevated and non-elevated position)
Ground Clearance <sup>2</sup>	1" at motor bracket
Front Wheels	6" flat free
Drive Wheels	10" flat free
Rear Wheels	6" flat free
Drivetrain	Two motor, Mid-Wheel 6 <sup>®</sup> drive
Braking System	Regenerative and electromechanical
Suspension Type	Active-Trac <sup>®</sup> (ATX) with Mid-Wheel 6 <sup>®</sup> Technology
High-back Seat Sizes	16"x16"-18", 18"x18"-20", 20"x20"-22"
Seat-to-Floor Range	21.25" – 33"
Per Charge Range <sup>2,4</sup>	Up to 16.8 miles (U1) Up to 19.2 miles (40 Ah)
Standard Electronics	40A, PG VR2
Component Weights <sup>6</sup>	Base: 252 lbs. w/ U1 batteries 267 lbs. w/ 40 Ah batteries U1 Batteries: 24.5 lbs. 40 Ah Batteries: 32.5 lbs.
Battery Requirements <sup>3,5</sup>	12 volt, U-1 (standard), 40 Ah (optional)
Battery Charger	5 amp, off-board
Warranty	Frame: 5-year limited Drivetrain: 13-month Electronics: 13-month Battery: 13-month

### HCPCS Code: K0899

POWER MOBILITY DEVICE, NOT CODED BY DME PDAC OR DOES NOT MEET CRITERIA

- Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition and tire condition. These specifications can be subject to a variance of (+/- 10%).
- Due to manufacturing tolerances and continual product improvement, this specification can be subject to a variance of (+ or - 3%).
- AGM or gel cell type required.
- Tested in accordance with ANSI/RESNA, WC Vol2, section 4 & ISO 7176-4 standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Test conducted at maximum weight capacity.
- Battery weight may vary based on manufacturer.
- Weight includes base, seat and batteries.

**PRIDE**  
MOBILITY  
LIVE YOUR BEST<sup>®</sup>

(US) 800-800-8586 • (Canada) 888-570-1113

pridemobility.com |    

HCPCS codes provided should not be considered as legal advice and do not guarantee reimbursement. DME providers are responsible for determining appropriate billing codes when submitting for insurance reimbursement.

The information contained herein is correct at time of publication; we reserve the right to alter specifications without prior notice.

NOTE: All specifications subject to change without notice. The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed and range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3%. Special order seats or non-standard dimensions will not be fold-flat and depth-adjustable.

JAZZYAIR2\_06-29-20

**Jazzy**  
air<sub>2</sub><sup>®</sup>

THE **FUTURE**  
OF **MOBILITY**



AMERICA'S FAVORITE POWER CHAIR

gojazzy.com |    

# EMPOWERING



Height and angle adjustable headrest

High-back or solid seat pan, reclining sport seat. Seat features 2" of depth adjustment. Choice of 16"x16"-18", 18"x18"-20", 20"x20"-22"

Newly re-designed lift features 12" of power adjustable seat height in 11 seconds



Convenient "Air" switch on armrest activates elevation

Front under seat storage keeps your valuables accessible



Active-Trac® Suspension with Mid-Wheel 6® Technology allows for maximum stability in the elevated and non-elevated positions.



Bright LED marker lights offers visibility in all environments.

## Elevate your life to the next level!


At the touch of a switch, the Jazzy Air® 2 elevates 12" in just 11 seconds. With a sleek, more modern look, this power chair safely drives at speeds up to 4 mph while elevated, allowing users to enjoy face-to-face social engagement. With Mid-Wheel 6® Drive and Active-Trac® Suspension, the Jazzy Air 2 offers enhanced performance and stability both indoors and out.

- Angle, width, height and depth adjustable armrests
- Integrated foot platform with rubber insert and 2" of adjustment for comfort in both the elevated and non-elevated positions
- Single-sided caster forks for enhanced maneuverability

### Available Colors



### Matte Colors

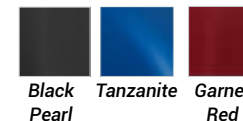

### Accessories

- Cane holder
- Cell phone holder
- Cup holder
- Flag kit
- 50", 60", 70" Lap belt
- Oxygen tank holder
- Rear basket
- Saddlebag
- Walker holder
- Weather cover
- XLR USB Charger

### Options

- 40-amp AGM batteries
- Swing-away joystick
- Stealth® cushions

### Military Patches

