

S I G N A T U R E S E R I E S



ALL NEW
Custom
Bucket Seat
& USB Port
Standard!

PR405-MLA Orion
with Twilight Positioning
shown in Anchor

The ONLY
3-position chair
that can lift your
feet above your
heart!



TWILIGHT



* upcharge applies

Orion
PR405-MLA

Introducing the world's first three-position power lift recliner with our patented Twilight Technology! Our dual-motor system comes standard and allows for a deeper recline to achieve rejuvenating positions no other three-position chair on the market offers. Rest, relax and rejuvenate in the Orion with Twilight for a weightless experience!



PR405-MLA

Model Number	PR405-MLA
Size	Medium/Large
Weight Capacity	375 LBS
Positioning	4-Way
Twilight Technology	YES
Chaise Seat	YES
Pocketed Coil Seat	NO
Back Type	Seam
Removable Back Option	NO
Head Rest Cover	None
Arm Rest Covers	None
Overall Width	40.5"
Overall Height	46.5"
Width Between Arms	23"
Floor to Top of Seat	21"
Seat Depth	21"
Seat to Top of Back	30.5"
Seat to Top of Arm	5"
Seat Angle Raised	17°
Floor to Top of Foot Rest Reclined	25"
Distance Required from Wall Reclined	21"
Extension reclined from Back of Seat to end of Foot Rest	43"
Life Time Warranty, Chair Frame, Lift Frame, Recline Mechanism	YES
Warranty on Electrical Parts	3 Yrs.
Pro Rated Years 4 Thru 7	YES
Warranty on Mechanical Labor	3 Yrs.
Warranty on Electrical Labor	1 Yr.
Customizable*	NO
Hand Control Routing	Inside Arm
Extra Pocket Option	Standard
Heat & Massage Option	NO
Power Pillow Option	NO
Power Adjustable Lumbar Option	NO
AutoFlex Option Available	NO
Foot Rest Extension Option	NO
Straight Lift Option	NO
USB Charging Port	YES
Place of Origin	USA



TWILIGHT




Motor #1

Controls the movement of both the backrest and footrest.

Motor #2

Activates Twilight technology, a tilting movement that provides a weightless experience.

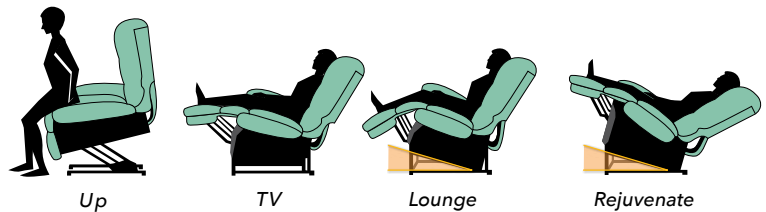

* See Customize Chart for Dimension limits
 Sizes may vary depending on fabric, filling material, upholstery and carpet thickness. All measurements taken on concrete floor using levels and metal rulers.



Smart Tek

In the event of a malfunction, our patented SmartTek™ diagnostic system enables you to get your lift chair up and running quickly.



Battery Backup System

Batteries not included. Battery Backup System is designed to bring the chair from the reclined to the seated position in the event of a power outage.

Literature is current at time of printing. Golden Technologies reserves the right to make changes to the product or literature at any time.



Complies with
 UL 73 and
 CSA C22.2 No. 68

